


carilex[®]

CARICHAIR

Advanced Cushion System Enables Mobility with Active Care

CariChair is a semi-automatic mobile cushion system with a battery backup, providing patients with uninterrupted care whenever and wherever they need it. With the leading-edge IPS technology, CariChair responds to patient movements and automatically adjusts internal cushion pressure.

CariChair is designed to be compatible with various mobile devices, including geriatric chairs, recliners, and wheelchairs, making it the ideal mobile system for comprehensive wound care therapy.

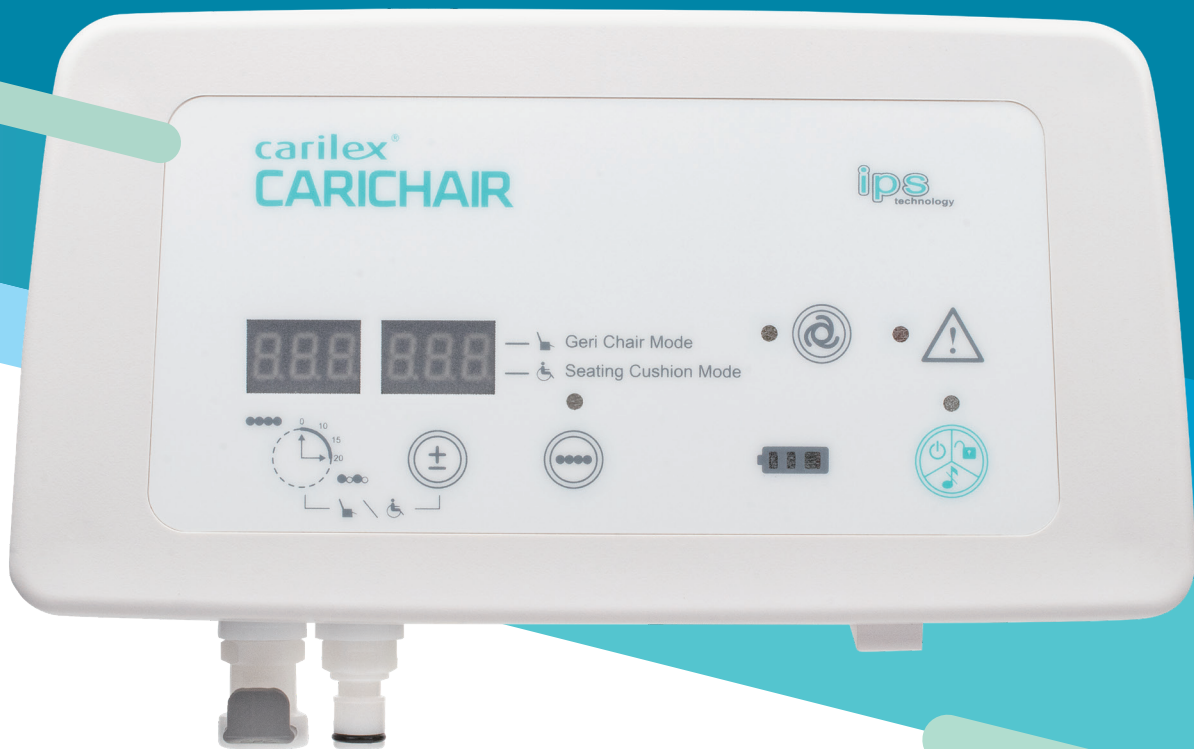

Intelligent Pressure Sensing (IPS) Technology

Responds to patient movements on mattress by adjusting internal mattress pressure automatically



Dual Therapy Modes

Provides both Dynamic and Static functions to meet individual needs



Comfort Setting

Provides optimal internal pressure according to each patient's sensation



Battery Backup up to 8 hours

Automatically activates battery backup in the event of main power failure and allows treatment to continue during patient transportation and examination

The CariChair Cushion

CariChair fits seamlessly with recliners and mobile devices, offering pressure relief and redistribution for patients with sensitive tissues.



Whole Body Therapy

The CariChair Cushion provides head-to-toe protection for patients even in sitting position.



Extra Support at the Sacrum Area

Specially designed deep air cells, for excellent immersion and pressure relief.



Quick Connect Coupling

Provides an easy connection between power unit and cushion.



Technical Data

Power Unit Model	REF : SR321
Power Input	100-240Vac 50-60Hz 1.5A
Power Consumption	23W Adaptor + Pump; 7W Pump only
Weight	Power Unit : 0.92Kg (2.0lb) Cushion Set : ● CariChair Cushion 2.1Kg (4.7lb) ● Mobicare Plus / Mobicare Cushion 1.5Kg (3.3lb)
Dimension	Power Unit : 10" x 6" x 2" (24.5 x 15.5 x 6.1cm) (W x H x D) Cushion Set : (W x L x H) ● CariChair Cushion, 71" x 20" x 6" (180 x 50 x 14 cm) ● Mobicare Plus Cushion, 18" x 18" x 4" (45 x 45 x 9 cm) ● Mobicare Cushion, 18" x 18" x 4" (45 x 45 x 9 cm)
Agency Approval	US NRTL - MET classified Mark_ANSI/AAMI ES60601-1
Operating Condition	● Temperature : 5°C to 40°C (41°F to 104°F) ● Humidity : 15% to 90%
Transportation & Storage Condition	● Temperature : -25°C to 70°C (-13°F to 158°F) ● Humidity : 0% to 90%
Weight Capacity	● CariChair PU+CariChair C-set: 150Kg (330lb) ● CariChair PU+Mobicare Plus set: 150Kg (330lb) ● CariChair PU+Mobicare C-set: 120Kg (260lb)

This Power Unit is not AP/APG protected

www.carilexmedical.com

Carilex Medical, Inc. **HO**

TEL : +886-3-328-7882
No.77, Keji 1st Rd., Guishan District,
Taoyuan City (33383), Taiwan

Carilex Medical USA

TEL : +1-866-551-4227
17068 Evergreen Place, Unit A, City of Industry,
CA 91745 USA

Carilex Medical Germany

TEL : +49-4121-476896
Fritz-Strassmann Str. 5, 25337 Elmshorn,
Germany

Carilex Medical United Kingdom

TEL : +44-1323-761900
5th Floor, 167-169 Great Portland St,
London W1W 5PF, United Kingdom